



WinForth



RF Therapy Platform Device!

BRIEF INTRODUCTION:

WinForth THERAPY DRAMATICALLY IMPROVES the RESULTS of your massage by reaching deep-seated conditions that hitherto have been deemed unreachable. The energy developed by your WinForth machine combines easy handling and intuitive functionality. You concentrate on your manual techniques and the energy compliments each of your expert movements, providing a treatment that is unparalleled in patient comfort and flexibility.

Provide your patients with a faster, more comfortable and pain free recovery.

ELIMINATES PAIN

WinForth energy inhibits pain immediately for a duration of more than 48 hours after treatment.

ACCELERATES HEALING

WinForth energy accelerates the body's natural self-repair system by promoting intra and extracellular exchanges in a lasting way. The benefits of your treatment continue after the session.

FREES UP MOVEMENT

WinForth energy relieves tension (contractures, oedema, fibrosis) by revascularising them locally. Your range of movement is immediately restored.

LOOK AESTHETIC

WinForth RF energy transform to thermal effect, speed up fat metabolism to reach body shaping results; stimulate collagen fibers to regrow and reconstruct, anti-aging, skin lifting and tightening.


HOW IT WORKS?

WinForth energy is a high frequency current oscillating between 300 kHz, 500kHz and 1 MHz. It accelerates "natural regeneration" of biological tissues. Its low intensity energy means it is completely non-invasive and 100% natural for the human body.

WinForth energy encourages intra and extra cellular exchanges and therefore faster tissue healing.

The natural resistance of tissues transforms WinForth energy into heat (diathermy). This localized increase in heat forces revascularization of restricted areas (fibrosis, contractures) and frees up joint movement.

WinForth THERAPY can act at different tissue depths. It has established 3 working frequencies to reach superficial pathologies, semi deep or deep conditions.

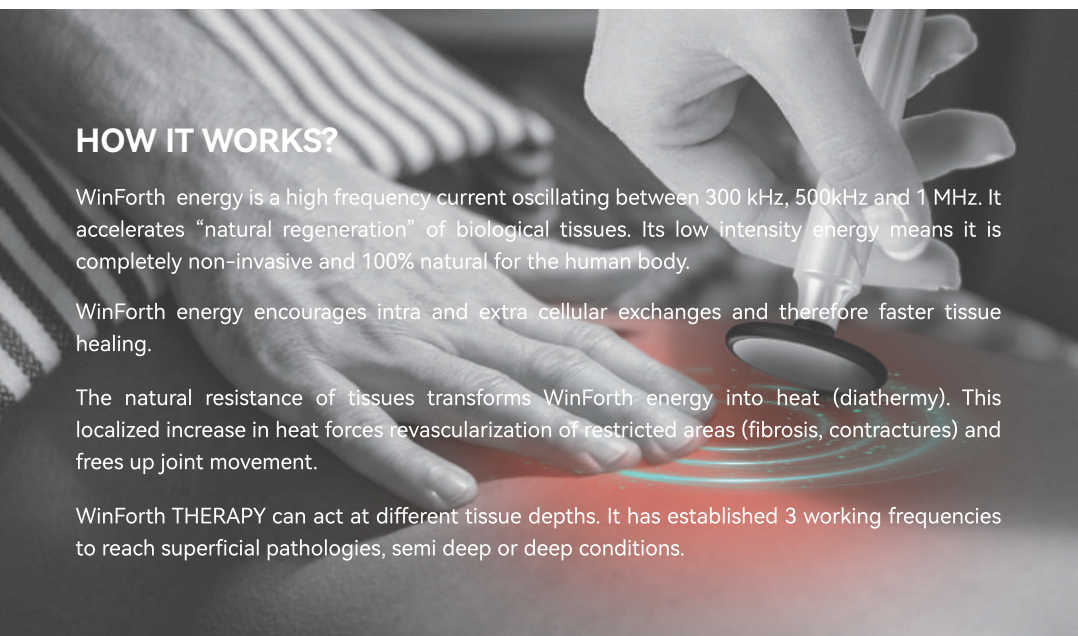



1 MHZ SUPERFICIAL ACTION

Skin layer heating

300 KHZ DEEP ACTION

HARD TISSUE = RET concentrates energy in all hard tissues such as joints, tendons and bones. The action is polyarticular and deep. LIGAMENT, JOINT, TENDON, PERIOSTEUM, BONES.

By wristband, your hand replaces the mobile electrode. The energy passes into the patient's body through your hand.

Completely non-invasive, this technique offers an exceptional experience where you can feel the effects of the energy directly, like a contracture easing within seconds in your hands.

SWAP MODE:

SUPERFICIAL ACTION/SEMI DEEP ACTION


1 MHz and 500 KHz RF energy Alternate output

SEMI DEEP ACTION/ DEEP ACTION

500 KHz and 300 KHz RF energy Alternate output

500 KHZ SEMI DEEP ACTION

SOFT TISSUE = CET concentrates energy in the muscles and other fluid tissues. The action is local, superficial or semi deep (DEEPCET) MUSCLE, FASCIA, APONEVROSIS, VISCERAL, LYMPHATIC AND VASCULAR.

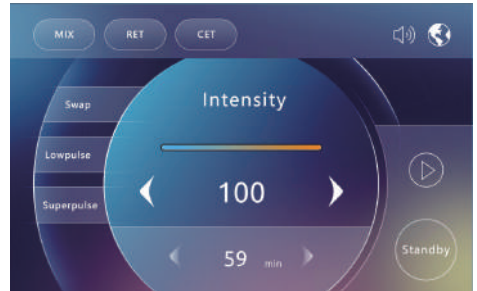

SUPER PULSE:


MAXIMUM POWER TO RELIEVE TENSION

Deeper thermal action.

LOW PULSE:

INFLAMMATION AND LYMPH NODES

More comfortable treatment.

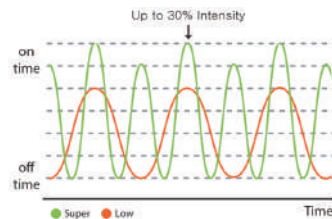

SPECIFICATION:

Power: MIX 35W; CET 165W; RET 165W

Frequency: 1 Mhz, 500 Khz , 300 Khz

Function: MIX; RET; CET

Modes: Swap; Low pulse; Super pulse

Accessories: 3 RET tips(30mm,60mm,80mm); 5 CET tips(20mm, 30mm, 40mm, 60mm, 80mm); MIX tip 39mm ; Negative plate; Negative pad; Wristband

Power Supply Input: 100-240VAC, 50-60Hz

Weight: 14Kg (Mainframe); 18Kg (Base)

Size: 490mm*430mm*260mm(Mainframe); 490mm*430mm*1020mm(Base)

APPLICATION:

Muscular, Joint And Lymph Problems

Pain Relief And Anti-Inflammatory; Assisted rehabilitation training

Traumatology

Rheumatology

Sports Physiotherapy; Enhanced sports training

Postpartum care: postpartum lochia discharge, mastitis, waist pain relief etc.

Fat metabolism, anti-aging, reducing cellulite, body shaping, skin tightening and lifting.



Shenzhen LeafLife Technology Co., Ltd.

National High-Tech Enterprise

ISO13485 Medical system certificate approved.

Innovative Benchmark Enterprise in the Medical Beauty Application Equipment Industry



Address:

4F,Bldg.C,JMD Industrial Park, No.39 Qingfeng Blvd., Baolong Industrial Area, Longgang Dist., Shenzhen, China

Email: info@leaflife.cn / info@leaflifetech.com

Website: www.leaflifetech.com

Tel.: 0755-27216609; 18682388756

